

**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

Data Requirement: PMRA Data Code:  
EPA DP Barcode: D354164  
OECD Data Point:  
EPA Guideline: Non-Guideline

**Test material:**

Common name: Metrafenone

Chemical name:

IUPAC name: (3-Bromo-6-methoxy-2-methylphenyl)(2,3,4-trimethoxy-6-methylphenyl)-methanone.

Benzophenone, 3'-bromo-2,3,4,6'-tetramethoxy-2',6-methyl-

3'-Bromo-2,3,4,6'-tetramethoxy-2',6-dimethylbenzophenone.

5-Bromo-6,6'-dimethyl-2,2',3',4'-tetramethoxy benzophenone.

CAS name: (3-Bromo-6-methoxy-2-methylphenyl)-(2,3,4-trimethoxy-6-methylphenyl)-methanone.

3'-Bromo-2,3,4,6'-tetramethoxy-2'6-methyl-methanone.

Methanone, (3-bromo-6-methoxy-2-methylphenyl)(2,3,4-trimethoxy-6-methylphenyl)-.

CAS No.: 220899-03-6.

Synonyms: AC 375839, AC375839, CL 375839, CL375839, BAS 560 F, BAS 560F, BAS 560 00 F, Benzophenone, ROI-24, Reference Substance 1, BAS560, Reg. No. 4037710, Aladin request number 27274.

Smiles string: C1=CC(Br)=C(C)C(C(=O)C2=C(C)C=C(O)C(O)=C2C)=C1C\_O\_CC\_C\_O\_CC\_CC\_C\_C\_C (EPI Suite, v3.12 SMILES string from ISIS .MOL file).

Primary Reviewer: Lynne Binari  
Cambridge Environmental

Signature:  
Date: 09/25/2008

Secondary Reviewer: Kathleen Ferguson  
Cambridge Environmental

Signature:  
Date: 09/25/2008

QC/QA Manager: Joan Gaidos  
Cambridge Environmental

Signature:  
Date: 09/25/2008

Final Reviewer: Lucy Shanaman  
EPA

Signature: *Lucy Shanaman*  
Date: 01/02/09

Company Code:

Active Code:

Use Site Category:

EPA PC Code: 000325.

**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

**CITATION:** Xing, J. 2002. BAS 560 F (CL 375839): laboratory validation of LC/MS determinative and LC/MS/MS confirmatory method M 3503 for the determination of BAS 560 F and CL 375816 residues in drinking and surface water. Unpublished study performed by BASF-Agro Research, Princeton, New Jersey; sponsored and submitted by BASF Corporation, Research Triangle Park, North Carolina. BASF Study No.: 92533, Report No.: RES 02-009 and Registration Document No.: 2002/7005048. Experimental start date January, 30, 2002, and completion date February 13, 2002 (p. 5). Final report issued March 20, 2002.

**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

**EXECUTIVE SUMMARY**

A laboratory validation of a method (BASF Agro Research Method M 3503) used to detect and quantify (3-bromo-6-methoxy-2-methylphenyl)(2,3,4-trimethoxy-6-methylphenyl)-methanone (metrafenone, BAS 560 F) and its transformation product 2-methyl-3-bromo-6-methoxybenzoic acid (CL 375816) in drinking and surface water was conducted to support registration of metrafenone (p. 5; Appendix B, pp. 56-58). This method validation was conducted in compliance with USEPA GLP Standards 40 CFR, Part 160 (p. 3). Drinking and surface water samples (see Table 1 below for characterization) were obtained on site at the test facility, with no additional details regarding collection provided (p. 8). Surface water was filtered (GF/B Grade, glass microfiber filter paper, double-layered) prior to use (Appendix B, pp. 61, 68). Aliquots (100 mL) of the waters were fortified with metrafenone (purity 99.7%, Lot No.: AC11957-109B) and CL 375816 (purity 95%, Lot No.: AC11238-127), in methanol, at 0.05 and 0.50 µg a.i./L (ppb, five treated samples per application rate and water type; pp. 8-9; Appendix B, pp. 62, 68). Metrafenone and CL 375816 were applied as a mixed fortification solution; the application solution volume used was not specified, but reported as ≤1.0 mL (Appendix B, pp. 62, 69). Ammonium hydroxide (1 mL) was added to the fortified water sample, then the sample was partitioned once with methylene chloride (100 mL; Appendix B, p. 69). The organic phase, containing metrafenone, was evaporated to dryness by rotary evaporation (*ca.* 35°C), with the resulting residues reconstituted in methanol (1.0 mL), using sonication, and diluted 1:1 (v:v) with water (1.0 mL) for LC/MS analysis (Appendix B, p. 69). The remaining aqueous phase was treated with acetic acid (4 mL), then partitioned twice with methylene chloride (100 mL x 2; Appendix B, p. 69). The organic phases, containing CL 375816, were combined and evaporated to dryness/near dryness as described above (Appendix B, p. 69). Resulting CL 375816 residues were reconstituted in methanol (*ca.* 5 mL) and taken to dryness; this step was repeated to completely remove trace acetic acid, then the residues were reconstituted in methanol (2 mL), using sonication, and diluted 1:1 (v:v) with water (2 mL; Appendix B, pp. 69-70). The CL 375816 sample was applied to a solid-phase extraction (SPE) cartridge (Isolute SAX, anion-exchange, 1 g/6 mL) preconditioned with methanol (*ca.* 5 mL), followed by methanol:water (1:1, v:v, *ca.* 5 mL; Appendix B, pp. 60, 70). The sample flask was rinsed with methanol:water (1:1, v:v, 2 mL) and the rinsate also applied to the cartridge. The loaded SPE cartridge was rinsed once with methanol (*ca.* 5 mL), then CL 375816 residues were eluted with methanol:formic acid (99:1, v:v, Solution C, *ca.* 5 mL; Appendix B, pp. 59, 70). The CL 375816 eluate was taken to dryness under nitrogen using a Mini-Vap (*ca.* 40°C) and the resulting residues reconstituted in water:methanol:acetic acid (69:30:1, v:v, Solution B, 1.0 mL) for LC/MS analysis (Appendix B, pp. 59, 70-71).

Metrafenone. Samples were analyzed by reverse-phase LC/MS under the following conditions: TosohHaas (currently Tosoh Bioscience LLC) TSK-GEL Super-ODS LC column (4.6 x 100 mm, 2 µm), guard filter and holder (TosohHaas No.: 18207 and 18206, respectively), gradient mobile phase combining (A) 0.05% aqueous acetic acid and (B) 10mM ammonium acetate in methanol [percent A:B at 0.0 min. 70:30 (v:v), 6.0-8.0 min. 0:100, 9.0-10.0 min. 70:30], column temperature ambient, injection volume 100 µL, flow rate 1.0 mL/minute, retention time 6.2-6.8 minutes, Finnigan MAT LCQ Deca MS with atmospheric pressure ionization interface, atmospheric pressure chemical ionization (APCI) positive ion mode, capillary temperature 150°C, vaporizer temperature 400°C, sheath gas flow 80 mL/minute, auxiliary gas flow 20 mL/minute, source voltage 6 kV, source current 5 µA, capillary voltage 41 V, MS full scan mode at m/z 200<sup>+</sup>-500<sup>+</sup> and SIM scan mode m/z

**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

411<sup>+</sup>/IsoW = 1.0 (pp. 9-11; Appendix B, pp. 64-65). As a second confirmatory method, samples were analyzed by LC/MS/MS as described above with the following modifications: MS/MS full scan mode m/z 150<sup>+</sup>-450<sup>+</sup>, SIM scan mode m/z 411<sup>+</sup>/IsoW = 1.5 (pp. 11-12).

CL 375816. Samples were analyzed by reverse-phase LC/MS as described above for metrafenone with the following modifications: gradient mobile phase [percent A:B at 0.0 min. 70:30 (v:v), 5.0 min. 45:55, 6.0-9.0 min. 0:100, 10.0-12.0 min. 70:30], flow rate 0.75 mL/minute, retention time 4.7-5.2 minutes, APCI negative ion mode, source voltage 5 kV, source current 80  $\mu$ A, capillary voltage -16 V, MS full scan mode m/z 180<sup>-</sup>-250<sup>-</sup>, SIM scan mode m/z 243<sup>-</sup>/IsoW = 1.0 (pp. 9-11; Appendix B, pp. 65-66). Samples were analyzed by LC/MS/MS with the following modification: MS/MS full scan mode m/z 180<sup>-</sup>-220<sup>-</sup> (pp. 12-13).

Metrafenone and CL 375816 extracted from the fortified samples were identified against reference standards (Figure 3, p. 31; Figures 6-8, pp. 34-36; Figures 11-13, pp. 39-41; Figures 15-16, pp. 43-44; Figure 18, p. 46). The limits of detection (LOD) and quantitation (LOQ) were 0.010 and 0.050  $\mu$ g a.i./L (ppb), respectively, for both compounds (Appendix B, p. 57).

Overall recoveries from water fortified at 0.05 and 0.5  $\mu$ g a.i./L averaged (n = 20) 101.0  $\pm$  3.4% (range 93.8-109.4%) of the applied for metrafenone and 83.0  $\pm$  10.0% (66.6-104.4%) for CL 375816 (DER Attachment 2). For each water type, recoveries for metrafenone and CL 375816 averaged (n = 10) 101.2  $\pm$  2.7% (98.4-105.2%) and 81.0  $\pm$  10.4% (68.8-104.4%), respectively, for the drinking water, and 100.8  $\pm$  4.0% (93.8-109.4%) and 85.0  $\pm$  9.1% (66.6-96.2%), respectively, for the surface water. For both water types, recoveries were comparable at the two fortification levels for metrafenone; however, average recoveries of CL 375816 were lower at the 0.05  $\mu$ g a.i./L fortification level (drinking water 76.9  $\pm$  6.9%, surface water 77.8  $\pm$  7.2%) as compared to the 0.5  $\mu$ g a.i./L level (drinking water 85.1  $\pm$  11.7%, surface water 92.2  $\pm$  2.8%). Detector responses were linear for metrafenone ( $r^2 = 0.9996$ ) over a range of 1.25-5.0  $\mu$ g/L and for CL 375816 ( $r^2 = 0.9998$ ) over a range of 2.5-10.0  $\mu$ g/L (p. 13; Figures 1-2, pp. 29-30). Chromatograms of unfortified drinking water detected apparent residues of metrafenone and CL 375816 at <0.00194  $\mu$ g a.i./L and <0.00445  $\mu$ g a.i./L, respectively, and at <0.00201  $\mu$ g a.i./L and <0.00690  $\mu$ g a.i./L, respectively, in unfortified surface water (pp. 5, 21-28; Figure 5, p. 33; Figure 10, p. 38; Figure 14, p. 42; Figure 17, p. 45). Minor interferences for both compounds were also detected in reagent blank samples (Figure 4, p. 32; Figure 9, p. 37); quantitative results were not provided.

The study author reported that processed water samples were determined to be stable for at least 7 days prior to analysis under refrigerator storage (ca. 4  $\pm$  2°C); however, no supporting data were provided (p. 6).

**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

Table 1: Properties of the waters.<sup>1</sup>

Property	Control Drinking Water	Control Surface Water
Sample ID No.	AC12584-10	AC12584-12
pH	6.1	5.6
Dissolved organic carbon (mg/L)	Not determined.	
Total organic carbon (mg/L)	0.8	5.9
Hardness (mg CaCO <sub>3</sub> /L)	102	47
Electrical conductivity (mmhos/cm)	0.42	0.19
Sodium Adsorption Ratio	0.94	0.69
Total Dissolved Solids (mg/L)	196	132
Turbidity (NTU)	1.32	4.92

Data obtained from p.8, Appendix C, pp. 98-99 of the study report.

<sup>1</sup> Both water samples were obtained on site at BASF Agro Research, Princeton, New Jersey (p. 8); no additional details regarding collection were reported.

**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

---

**Attachment 1: Structures of the Test Material**

**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

**Metrafenone [AC 375839, AC375839, CL 375839, CL375839, BAS 560 F, BAS 560F, BAS 560 00 F, Benzophenone, ROI-24, Reference Substance 1, BAS560, Reg. No. 4037710, Aladin request number 27274]**

**IUPAC Name:** (3-Bromo-6-methoxy-2-methylphenyl)(2,3,4-trimethoxy-6-methylphenyl)-methanone.

Benzophenone, 3'-bromo-2,3,4,6'-tetramethoxy-2',6-methyl-  
3'-Bromo-2,3,4,6'-tetramethoxy-2',6-dimethylbenzophenone.  
5-Bromo-6,6'-dimethyl-2,2',3',4'-tetramethoxy benzophenone.

**CAS Name:** (3-Bromo-6-methoxy-2-methylphenyl)-(2,3,4-trimethoxy-6-methylphenyl)-methanone.


3'-Bromo-2,3,4,6'-tetramethoxy-2'6-methyl-methanone.

Methanone, (3-bromo-6-methoxy-2-methylphenyl)(2,3,4-trimethoxy-6-methylphenyl)-.

**CAS Number:** 220899-03-6

**SMILES String:** CC1=CC(Br)=C(C)C(C(=O)C2=C(C)C=C(OC)C(OC)=C2C)=C1COCOC  
CC1=CC(C)=C(OC)C(OC)=C1C(=O)C2=CC(Br)=C(C)C(OC)=C2 (EPI Suite, v3.12 SMILES string from ISIS.MOL file).

**Empirical formula:** C<sub>19</sub>H<sub>21</sub>BrO<sub>5</sub>      **Molecular formula:** C<sub>19</sub>H<sub>21</sub>BrO<sub>5</sub>

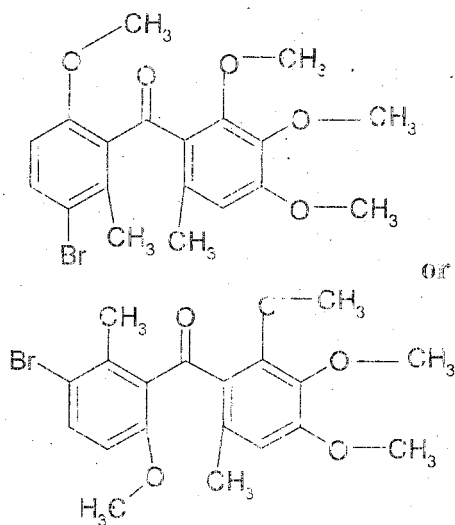

**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

**AC 375816 [CL 375816, ROI-7, Reference Substance 3, AQ, Reg. No. 375816, Reg. No. 4074484, Aladin request number 27100]**

**IUPAC Name:** Not reported.


**CAS Name:** 2-Methyl-3-bromo-6-methoxybenzoic acid.  
3-Bromo-6-methoxy-2-methylbenzoic acid.  
Benzoic acid, 3-bromo-6-methoxy-2-methyl-.  
Benzoic acid, (2-methyl-3-bromo-6-methoxy)-.

**CAS Number:** Not reported.


**SMILES String:** CC(=C(Br)=C(C)C(=O)OC)C (EPI Suite, v3.12 SMILES string from ISIS .MOL file).

**Empirical formula:** C<sub>9</sub>H<sub>9</sub>BrO<sub>3</sub>

**Molecular formula:** C<sub>9</sub>H<sub>9</sub>BrO<sub>3</sub>

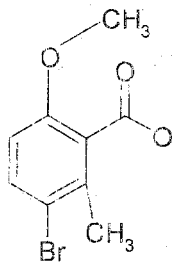

or

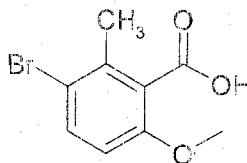



**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

---

**Identified Compounds**

**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

**Metrafenone [AC 375839, AC375839, CL 375839, CL375839, BAS 560 F, BAS 560F, BAS 560 00 F, Benzophenone, ROI-24, Reference Substance 1, BAS560, Reg. No. 4037710, Aladin request number 27274]**

**IUPAC Name:** (3-Bromo-6-methoxy-2-methylphenyl)(2,3,4-trimethoxy-6-methylphenyl)-methanone.

Benzophenone, 3'-bromo-2,3,4,6'-tetramethoxy-2',6-methyl-

3'-Bromo-2,3,4,6'-tetramethoxy-2',6-dimethylbenzophenone.

5-Bromo-6,6'-dimethyl-2,2',3',4'-tetramethoxy benzophenone.

**CAS Name:** (3-Bromo-6-methoxy-2-methylphenyl)-(2,3,4-trimethoxy-6-methylphenyl)-methanone.


3'-Bromo-2,3,4,6'-tetramethoxy-2',6-methyl-methanone.

Methanone, (3-bromo-6-methoxy-2-methylphenyl)(2,3,4-trimethoxy-6-methylphenyl)-.

**CAS Number:** 220899-03-6

**SMILES String:** C1=CC(Br)=C(C)C(C(=O)C2=C(C)C=C(O)C(O)=C2C)=C1C\_O\_CC\_C\_O\_CC\_CC\_C\_C\_C (EPI Suite, v3.12 SMILES string from ISIS .MOL file).

**Empirical formula:** C<sub>19</sub>H<sub>21</sub>BrO<sub>5</sub>      **Molecular formula:** C<sub>19</sub>H<sub>21</sub>BrO<sub>5</sub>

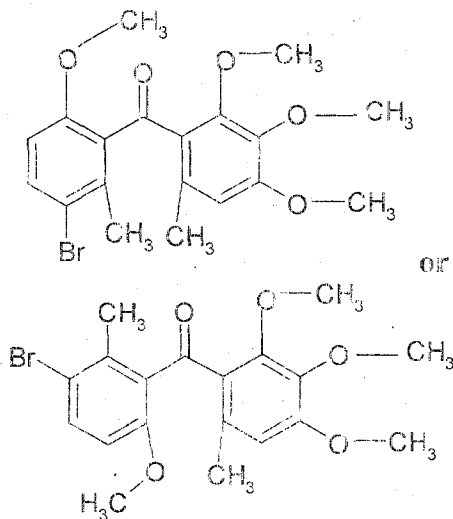

**Data Evaluation Record on method validation for metrafenone (BAS 560 F) and its product CL 375816 in water**

PMRA Submission Number {.....}

EPA MRID Number 47267436

**AC 375816 [CL 375816, ROI-7, Reference Substance 3, AQ, Reg. No. 375816, Reg. No. 4074484, Aladin request number 27100]**

**IUPAC Name:** Not reported.


**CAS Name:** 2-Methyl-3-bromo-6-methoxybenzoic acid.  
3-Bromo-6-methoxy-2-methylbenzoic acid.  
Benzoic acid, 3-bromo-6-methoxy-2-methyl-.  
Benzoic acid, (2-methyl-3-bromo-6-methoxy)-.

**CAS Number:** Not reported.

**SMILES String:** CC(=C(Br)C(C)C(=O)OC)C (EPI Suite, v3.12 SMILES string from ISIS .MOL file).

**Empirical formula:** C<sub>9</sub>H<sub>9</sub>BrO<sub>3</sub>

**Molecular formula:** C<sub>9</sub>H<sub>9</sub>BrO<sub>3</sub>

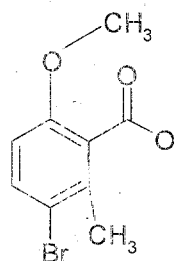

or

